

---

# Advanced Normalization Tools (ANTs)

*Release 1.5*

Brian B. Avants, Nick Tustison and Gang Song

March 20, 2011

Penn Image Computing And Science Laboratory  
University of Pennsylvania

## Abstract

We provide examples and highlights of Advanced Normalization Tools (ANTs) that address practical problems in real data.

## Contents

<b>1</b>	<b>Introduction</b>	<b>2</b>
1.1	The ANTS Executable	2
1.2	The WarpImageMultiTransform Executable	3
1.3	I/O Data Formats In ANTS	3
	Image Volumes	3
	Affine Transformation File	4
	Deformation Field File	4
	Labeled Point Sets	4
1.4	The ImageMath Executable	5
1.5	ANTS Scripts	5
<b>2</b>	<b>ImageRegistration with ANTS</b>	<b>6</b>
2.1	Use Your Header!	6
2.2	ANTS Transformation Models	6
	Affine and Rigid Registration	7
	Deformable Registration	9
2.3	ANTS Similarity Terms and Other Functions	13
2.4	Choosing a Metric	16
2.5	Diffusion Tensor Normalization	19
2.6	Multivariate Normalization with ANTS	19
2.7	Optimal Template Construction with ANTS	20
	A Concrete Template Construction Example, Step-by-Step	20
2.8	Notes on Large Deformation Mapping	21
<b>3</b>	<b>Image Segmentation and Labeling</b>	<b>24</b>
3.1	Basic Segmentation	24

---

3.2	Prior and Template-based Image Segmentation . . . . .	24
3.3	Cortical Thickness . . . . .	25
3.4	User/Label-Guided Normalization for Hippocampus Mapping . . . . .	25
<b>4</b>	<b>Application to Studies</b>	<b>27</b>
4.1	Brain Mapping in the Presence of Lesions . . . . .	27
4.2	Statistical Mapping with ANTS: Morphometry, Function, Jacobian, Thickness . . . . .	27
4.3	Statistics with ANTS and <b>R</b> . . . . .	29
<b>5</b>	<b>Dependencies and Related Software</b>	<b>32</b>
5.1	Dependencies and Compilation . . . . .	32
5.2	Visualization and Quantification of ANTS Results . . . . .	32
5.3	Statistics Beyond ANTS . . . . .	32
5.4	Pipelining with ANTS . . . . .	32
<b>6</b>	<b>Annotated Bibliography</b>	<b>33</b>

---

## 1 Introduction

The ANTS framework provides open-source functionality for deformable normalization with large deformations, as shown in figure 1. Small deformation mappings and segmentation tools are also available. Independent evaluation of ANTS normalization tools, applied to “control” data, placed the toolkit as a top performer amongst 14 methods [23]. Developer evaluation showed stronger differences with other methodology in neurodegenerative neuroimaging data, where large deformation is required [10]. ANTS serves as both a base for further algorithm development and also as an application-oriented toolkit. ANTS enables diffeomorphic normalization with a variety of transformation models, optimal template construction, multiple types of diffeomorphisms, multivariate similarity metrics, diffusion tensor processing and warping, image segmentation with and without priors and measurement of cortical thickness from probabilistic segmentations. The normalization tools, alone, provide a near limitless range of functionality and allow the user to develop customized objective functions. Objective functions in ANTS are of the form:

$$\text{Deformation Cost} + \text{Data Terms},$$

and the command line reflects this balance of two terms. As mentioned above, the data term may combine multiple different measures of similarity that are optimized in parallel, for instance, image similarity and landmark terms. This document seeks to provide a practical overview of basic functionality and some of the common use cases that users seek. Additional information is available online – see [ANTS Homepage : www.picsl.upenn.edu/ANTS](http://picsl.upenn.edu/ANTS). For compilation details, see: [Compile and download: http://picsl.upenn.edu/ANTS/download.php](http://picsl.upenn.edu/ANTS/download.php). The most important ANTS programs are described below.<sup>1</sup>

### 1.1 The ANTS Executable

The ANTS program itself is the central program encapsulating normalization/registration functionality. Its main output is an affine transform file and a deformation field, potentially with inverse. Options to ANTS

---

<sup>1</sup>This document is a work in progress. Please check for updates with each release.

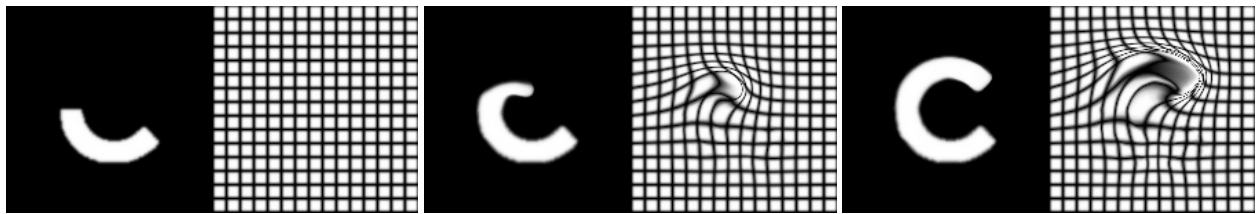


Figure 1: The original goal of ANTS was to develop public, open source large deformation image registration. This is a classic example showing the progress of deforming a half C to a full C along a geodesic diffeomorphism. The deforming grid accompanies each deformed image.

allows the user to navigate the similarity and transformation options that are available. ANTS takes multiple similarity and optimization criteria as options. The program is wrapped in `ants.sh` for normalization with “out of the box” parameters and in `buildtemplateparallel.sh` for computationally distributed optimal template construction. You can also call `ANTS --help` to get basic usage.

## 1.2 The WarpImageMultiTransform Executable

The `WarpImageMultiTransform` program applies ANTS mappings to images. One may apply an arbitrarily long series of transformations to an image through this program. Thus, it enables one to compose a series of affine and deformable mappings – and/or their inverses – such that one may avoid multiple interpolations of a single image. This program is used internally in ANTS optimization to avoid ever using more than one interpolation. The `use-NN` option applies nearest neighbor interpolation, otherwise linear interpolation is used. You can call `WarpImageMultiTransform--help` to get basic usage.

Also, please take a look at <http://picsl.upenn.edu/ANTS/ioants.php> and the section below on I/O.

## 1.3 I/O Data Formats In ANTS

ANTS support both 2D and 3D images. Since ANTS is implemented in ITK, it is able to read and write the popular data formats supported through ITK. There are four basic types of data for ANTS.

### Image Volumes

ANTS support both 2D and 3D images, including

- Nifti (.nii, .nii.gz)
- Analyze (.hdr + .img / .img.gz)
- MetaImage (.mha)
- Other formats through `itk::ImageFileWriter` / `itk::ImageFileReader` such as jpg, tiff, etc. See ITK documentation.

### Affine Transformation File

ANTS uses the `itk::TransformFileReader` / `itk::TransformFileWriter` to handle affine transformation files. These are text files recording the matrix and offset for the affine transform **FIXME NEED MORE DETAIL on definition and geometry of affine matrix**. The interpretation depends on the header but generally it breaks down like this:

```
#Insight Transform File V1.0
```

```
# Transform 0
```

```
Transform: MatrixOffsetTransformBase_double_3_3
```

```
Parameters: a b c d e f g h i j k l
```

```
FixedParameters: M N O
```

The meaning of these different terms:

The matrix component (rotation, scaling, shearing all composed together):

```
a b c
d e f
g h i
```

is followed by *j k l* — this is the translation.

*M N O* is the center of rotation .

You can use the `-initial-affine` and `-continue-affine` options to control how this is used in registration.

### Deformation Field File

ANTS writes the deformation field into a multi-component image (latest release). Thus, one refers to the warp file as `ABCD.nii.gz` which, in 3D, contains 3 components. If you want to see the details of how a warp is achieved, look at the code in `itkWarpImageMultiTransform.txx`. One issue to be aware of is that if you are directly indexing a vector component image in voxel coordinates, then you need to apply the image's direction matrix to the vector. An example of this is in `ANTSJacobian.cxx`, the `TransformVector` function.

### Labeled Point Sets

There are three formats supported

- VTK format: point + scalar label in a `.vtk` file. **FIXME NEED MORE DETAIL – what type of vtk data ? polydata?**.
- Labeled point sets are saved as an image volume with the intensity = 1,2,3 for label = 1,2,3. For example, figure 16 shows 2D images representing a 3-class labeled point set.
- Avants Format. In this simple format, the point sets are saved as a text file. The first and last line is always 0 0 0 0. Each line in between for a point as X Y Z L (L is the label.)

The names and formats of the output files are both specified by `-o` option to the ANTS executable. For example, for in 3D with `-o RegA2B.nii.gz` , the typical output files are: Affine transform file:

RegA2BAffine.txt Warping deformation field: RegA2BWarp.nii.gz. Inverse warping deformation field: RegA2BInverseWarp.nii.gz. To apply the affine transform and the deformation field to the moving image, use WarpImageMultiTransform. For example: to apply the transform to the moving image:

```
WarpImageMultiTransform 2 A.nii.gz ADeformed.nii.gz -R B.nii.gz RegA2BWarp.nii.gz
RegA2BAffine.txt
```

to apply the inverse transform to the fixed image:

```
WarpImageMultiTransform 2 B.nii.gz BDeformed.nii.gz -R A.nii.gz -i RegA2BAffine.txt
RegA2BInverseWarp.nii.gz
```

## 1.4 The ImageMath Executable

This is a multi-purpose program that has the following syntax: “ImageMath ImageDimension outputfilename Operation InputFileName parameters “. Most basic scalar image operations – and some tensor operations – may be performed with this program. Some operations output text files, some output images and some output only to the terminal. ImageMath allows one to multiply images together (m), to negate images (Neg), to take an image to a power (pow), to test the invertibility of transformations (InvId), to compute the fractional anisotropy of an image (TensorFA) and to perform many other feats of image processing strength. Like all other ANTS programs, one may call ImageMath from the command line to see all of its options. ImageMath is used heavily in ANTS scripts. You can call it, in compination with grep or less, to attempt to find some of its capabilities.

## 1.5 ANTS Scripts

The ANTS/Scripts directory contains (hopefully) user-friendly wrappings of ANTS tools that enable higher-level error checking and combinations of basic ANTS functions. These scripts are called as `bash antsscriptname.sh` and provide usage when called from the command line. For instance, try `bash ants.sh`. Additionally, if you run ants in tcsh shells, you will need to call it as follows:

```
~/Code/ANTSSourceForgeSVN/bin/ANTS 2 -m CC\[r16slice.nii,r64slice.nii,1,2\]
```

instead of

```
~/Code/ANTSSourceForgeSVN/bin/ANTS 2 -m CC[r16slice.nii,r64slice.nii,1,2]
```

because tcsh does not interpret the brackets as intended. Using the bash shell script interface also avoids this. The key change involves using the slash to force tcsh to interpet the bracket literally.

## 2 ImageRegistration with ANTS

There are two general applications of image registration. The first application is transforming labeled data from a template image space into an individual space. This strategy is important when appearance alone is not enough to locate a structure as, for example, in the case of hippocampus segmentation. The template provides a prediction of the hippocampus-amygdala boundary. The second application operates in the "inverse" direction of the first: instead of mapping template to individual, we map individual(s) to the template. Voxel-based population studies of either functional or structural variables depend on mapping to a template space. The common coordinate system enables a statistical evaluation of the likelihood of consistent activation across a group or, in other contexts, the differences in anatomy between two groups.

**Quick Start:** call `ants.sh` (from ANTS/Scripts/) to get usage and apply a normalization to some of your existing data. It is instructional to read the script, modify some of the parameters and re-run to witness the effect of your changes – image registration is an art as well as science. Many other "ready to go" scripts are available in ANTS/Scripts. The user must set his/her `ANTSPATH` environment variable – which points to the location of ANTS binaries – within these scripts or in their environment. Note: All ANTS programs provide usage when called from the command line. Most require the image dimension to be specified as the first parameter. E.g. `ImageMath ImageDimension` where `ImageDimension` is 2, 3 or 4.

The ANTS toolkit enables both types of mapping.

The main challenge in image and brain mapping is defining the way in which images/anatomy are compared. There are two components to the comparison. The shape transformation space defines the range of shape variation that will be allowed in the optimization. The appearance similarity space defines the statistical assumptions that determine when one image is considered to appear similar to another.

These two components interact in a weighted optimization within a multiple resolution gradient descent framework. Each component may use either "strict" or "flexible" assumptions about the shape or appearance similarity. The selection of these models should be done in a principled way that fits the problem at hand. No single choice is appropriate for all scenarios (see the "no free lunch" theorem).

Thus, ANTS enables many operating points from both the transformation and appearance domains such that users may make choices appropriate for their problems. The ANTS command-line syntax, shown in figure 2, reflects these operating points and the various components that interact in the optimization. ANTS may be used to navigate the transformation and similarity metric space.

### 2.1 Use Your Header!

ANTS uses the direction/orientation, spacing and origin in its definitions of the mapping between images. That is, ANTS uses the nifti standard definitions of an image space. So, the image orientation (direction), origin and spacing are important! The program 'PrintHeader' will give you information about your image's header.

### 2.2 ANTS Transformation Models

The ANTS toolkit provides a hierarchy of transformations with adjustable levels of complexity, regularization, degrees of freedom and behavior as optimizers. The simplest transformation model is the rigid and/or affine transform. The most complex – and most flexible – is a symmetric diffeomorphic transformation

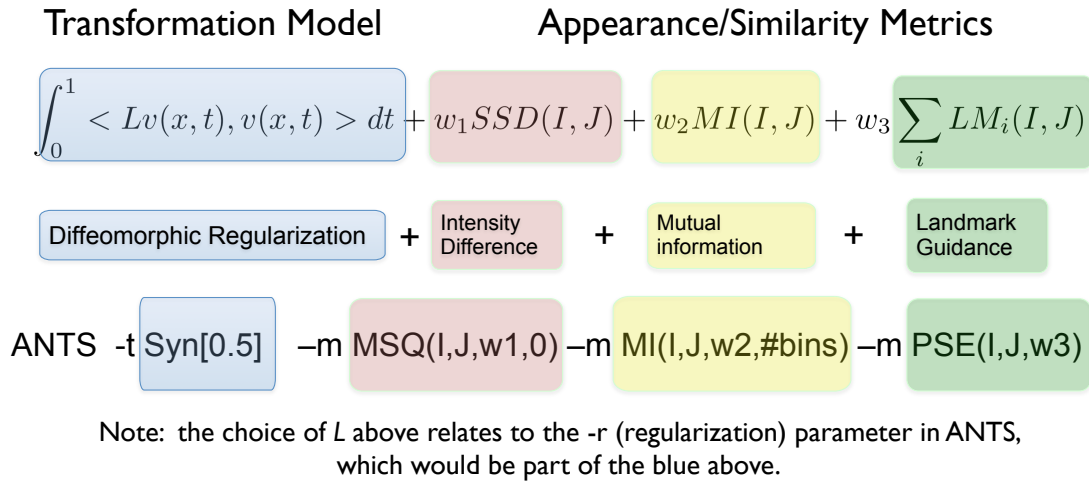


Figure 2: The relationship between the variational form that defines the optimization goals and the ANTS command line. The Syn option implements the method in [10] and evaluated in [23].

based on optimizing and integrating a time-varying velocity field. Computation time also increases with transformation model complexity, but most normalization needs may be met with under an hour of computation. We provide an overview of some of the ANTS models below and try to communicate intuition on what one gains/loses with each choice. An overview of available models and similarity terms is in Table 1.

### Affine and Rigid Registration

The most basic mapping between images is an affine mapping. The ANTS affine mapping syntax is shown in the affine mapping figure 3. The parameter breakdown of the ANTS call in the affine figure is:

- ANTS – the normalization program
- 2 – the expected image dimension
- -m MI[...,...,1,32] – the similarity metric “mutual information” with a 32-bin square joint histogram and weight of 1. Metric details will be discussed in a later section.
- -i 0 – no deformable iterations (affine only).
- -o ab – the output prefix. A specific extension may also be specified, e.g. -o ab.mhd .

The parameter breakdown of the WarpImageMultiTransform call in the affine figure is:

- WarpImageMultiTransform – to apply the mapping output from ANTS to an image.
- 2 – the expected image dimension.
- the image to be deformed is input next.
- the deformed image output file name is next.
- the transform itself is next – if is affine and preceded by “-i”, we apply the inverse affine map.

Category	Transformation, $\phi$	Similarity Measures	Brief Description
<b>Linear</b>	Rigid <sup>†</sup>	MI, MSQ	Rigid registration.
	Affine <sup>†</sup>	MI, MSQ	Affine registration.
<b>Elastic</b>	Deformable	CC, PR, MI, MSQ, PSE	Demons-like algorithm.
	DMFFD	CC, PR, MI, MSQ, PSE	FFD variant.
<b>Diffeo.</b>	Exponential <sup>†</sup>	CC, PR, MI, MSQ, PSE	$\min v(\mathbf{x})$
	Greedy SyN <sup>†</sup>	CC, PR, MI, MSQ, PSE	locally in time $\min v(\mathbf{x}, t)$
	Geodesic SyN <sup>†</sup>	CC, PR, MI, MSQ, PSE	$\min v(\mathbf{x}, t)$ over all time

Table 1: Transformations and a subset of the similarity metrics available in ANTS. Similarity metric acronyms: CC = fast cross correlation (the preferred metric), PR = cross correlation (deprecated), MSQ = mean squared difference, MI = mutual information, PSE = point-set expectation [25]. ANTS also provides the inverse of those transformations denoted by the ‘<sup>†</sup>’ symbol. The brief descriptions of the diffeomorphic algorithms contrast the way in which the velocity field is optimized and used to parameterize  $\phi$ , the mapping. All ANTS Diff algorithms generate  $\phi(\mathbf{x}, t)$  over  $t \in [0, 1]$  through gradient descent.

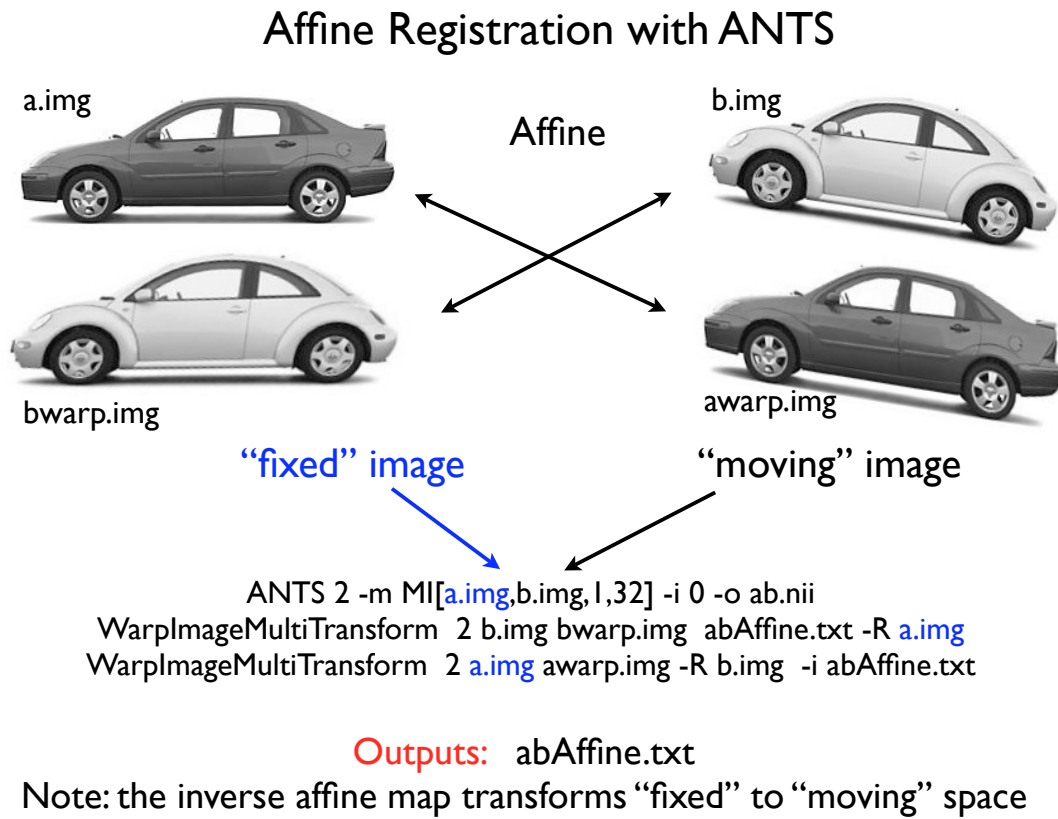


Figure 3: The anatomy of an ANTS optimization and application of the resulting warping.



- the ”-R ” option dictates the domain you want to warp into – usually the “fixed” image, unless using the inverse mapping, in which case one switches the role of fixed and moving images.

Note that -i 0 (in the ANTS call) means that no deformable mapping is performed. The call ”ANTS 2” means to treat the images as 2-dimensional. ”ANTS 3” would be called for 3D images. The ”-m ” term defines the similarity metric. Note that no more than one image interpolation is used during optimization. That is, ANTS always references the original image before applying a transformation or series of transformations.

Affine mapping may become nontrivial when the region of interest occupies only a small portion of the image data. In such cases, ANTS enables the user to define a mask to focus the optimization effort on the region of interest. Example code for (affine) registration with a mask:

```
ANTS 3 -m MI[AtlasHead.nii,crAnatomical.nii,1,32]
-o TEST -i 10x10x0 -r Gauss[1.5,0] -t Exp[0.5] -x mask.nii.gz
```

The difference between this command and a regular ANTS call is the -x option, which specifies the mask, defined in the template space (here, AtlasHead.nii). The mask option also affects the deformable optimization. Other affine registration options include `--number-of-affine-iterations 5000x5000x5000` which specifies that the affine registration uses a 3 level image pyramid with each level 5000 iterations at most. `--MI-option 48x6000` means to use Mutual Information as similarity metric with 48 bins and 6000 samples. `--affine-gradient-descent-option` defines the options for the gradient descent procedure; the 4 numbers are maximum step length, relaxation factor, minimum step length and translation scale.

**Rigid Registration:** ANTS will also perform rigid registration. To force the linear component of registration to be purely rigid (no scaling or shearing), then pass the parameter `--rigid-affine true` to the ANTS command line. For example,

```
ANTS 3 -m MI[AtlasHead.nii,crAnatomical.nii,1,32]
-o TEST -i 0 --rigid-affine true
```

The jacobian of the resulting affine mapping should be unity.

### Deformable Registration

Affine mapping is adequate when the difference between images involves only rotation, scaling or shearing. Other data requires more deformable mappings to capture shape differences and find a good alignment of image features. Figure 4 shows how deformable mapping may improve the correspondence of the deformed beetle to the ford image. Most ANTS-based applications use symmetric diffeomorphic normalization. However, ANTS also enables a simpler parameterization of a deformable mapping via a regularized vector space. We term these types of transformations as “elastic”. The original Demons algorithm provides a classic example of using a regularized vector space for nonlinear registration. Caveats are that a regularized vector space may not preserve the underlying topology and may also prove too inflexible to capture the shape changes one is after. Both of these shortcomings motivate the use of diffeomorphisms.

**Elastic/Vector Space Transformations.** If we assume no affine transformation, then an elastic transformation involves computing a mapping from image  $I(\mathbf{x})$  to image  $J(\mathbf{x})$  through a deformation field  $\mathbf{u}(\mathbf{x})$ . The deformation is defined in the physical space of the image and dictates the positional difference between corresponding features in the two images. Thus, if a feature defined at  $I(\mathbf{x})$  matches a feature in  $J$  at position  $\mathbf{y}$  then the deformation field at  $\mathbf{x}$  should give  $\mathbf{u}(\mathbf{x}) = \mathbf{y} - \mathbf{x}$ . Such a deformation field may be applied to



Figure 4: This example shows the degree to which the beetle (b.img) may be deformed to the ford (a.img) under different transformation models. Left to right increases the degrees of freedom in the mapping and thus the registration accuracy.

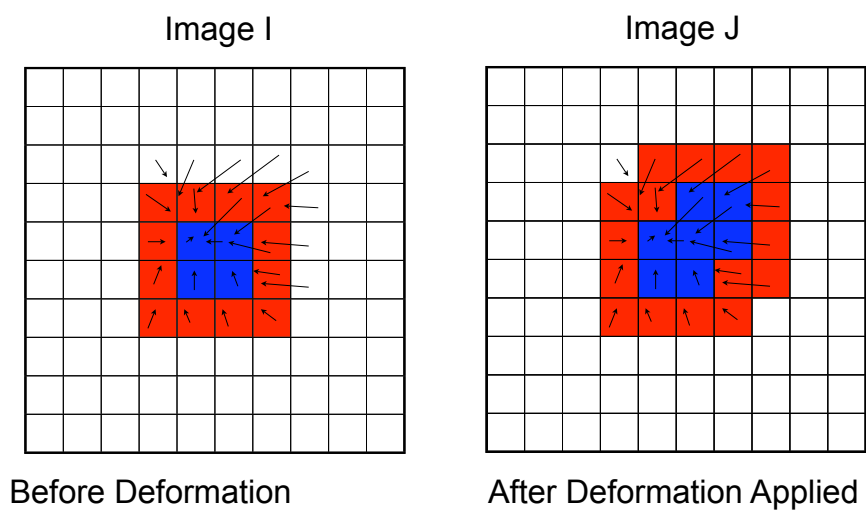


Figure 5: A deformation in a digital domain.

deform image  $J$  into image  $I$  by composing the mapping  $J_{deformed}(\mathbf{x}) = J(\mathbf{x} + \mathbf{u}(\mathbf{x}))$ . In a perfect world, then  $I(\mathbf{x}) = J_{deformed}(\mathbf{x})$ , though this is rarely the case. Figure 5 visualizes the deformation of an image under this standard model. Gradient descent optimization of an elastic mapping may be summarized (crudely) as:

$$\begin{aligned} &\text{Compute the similarity gradient: } \nabla E = \partial_{\mathbf{u}} \Pi(I, J(\mathbf{x} + \mathbf{u}(\mathbf{x}))). \\ &\text{Update the deformation field: } \mathbf{u}(\mathbf{x}) \leftarrow \mathbf{u}(\mathbf{x}) + \delta \nabla E. \\ &\text{Regularize the deformation field: } \mathbf{u}(\mathbf{x}) \leftarrow G_{\sigma} \star \mathbf{u}(\mathbf{x}), \end{aligned} \quad (1)$$

where  $\Pi$  is the similarity,  $\delta$  is a gradient step length and  $G_{\sigma}$  is a gaussian smoother. This optimization is captured in the following ANTS command:


```
ANTS 3 -m CC[AtlasHead.nii,crAnatomical.nii,1,4] -r Gauss[0,3]
-o ElasticTest.nii.gz -i 30x20x10 -t Elast[1.5]
```

Here, we use the CC metric (a cross-correlation implementation) with window radius 4, weight 1 and gradient step length 1.5. The optimization will be performed over three resolutions with a maximum of 30 iterations at the coarsest level, 20 at the next coarsest and 10 at the full resolution. We use a Gaussian regularizer with a sigma of 3 that operates only on the deformation field and not on the similarity gradient, as 0 is passed as the first parameter to Gauss. One may see the correspondence, yet again, between the ANTS call and the optimization scheme. The optimization will stop when either the energy cannot be smoothly minimized or the maximum number of iterations is reached. BSpline regularization is also available in ANTS – see the DMFFD section below.


**Warping and Invertibility.** ANTS uses physical space to define mappings. The origin etc is in the coordinates of the world in which the picture (e.g. MRI) was taken. Then, warping between physical coordinates is relatively easy. Differences in bounding boxes etc present no problem – assuming you avoid inconsistent headers i.e. origins, directions, data orientation. One may use PrintHeader to check the data and run simple, fast tests (few iterations) to perform sanity checks before running through loads of data. An ANTS deformation consists of a standard naming prefix plus a standard naming extension. We usually assume nii but other file types may be used. The standard naming extension is Warpxvec.nii.gz for the x component of the deformation with similar naming for y and z. Each component is stored in a separate scalar image. The value of a voxel of a deformation/warp component image is the physical space displacement from that voxel. Note that the inverse mapping is stored in InverseWarpxvec.nii.gz and provides the mapping in the opposite direction. Note that an inverse – in the digital domain – is only approximate as shown in figure 6. A few important notes follow. (1) Deformation directionality: Warps/deformations applied to images occur in the opposite direction from warps/deformations applied to points. That is, if we map image B to A using

```
ANTS 2 -m CC[A,B,1,2] -o OUTPUT
```

then we get a deformation called OUTPUTWarp+extensions and an affine transform called OUTPUTAffine.txt. This composite mapping - when applied to B - will transform B into the space of A. However, if we have points defined in B that we want to map to A, we have to use OUTPUTInverseWarp and the inverse of OUTPUTAffine.txt. This is because the domain and range of the map's definition need to be dense in a way that is appropriate for the data type. Figure 6 illustrates the concept. (2) Older Image Formats: older image formats (e.g. Analyze) may not have proper origin/offset/direction support! In these cases, we recommend converting to nii and verifying that data overlays properly. (3) Transform Composition: Composition of transforms may be achieved with ComposeMultiTransform. (4) Warping with inverses and concatenations – viable when using diffeomorphisms – are described in figure 7.

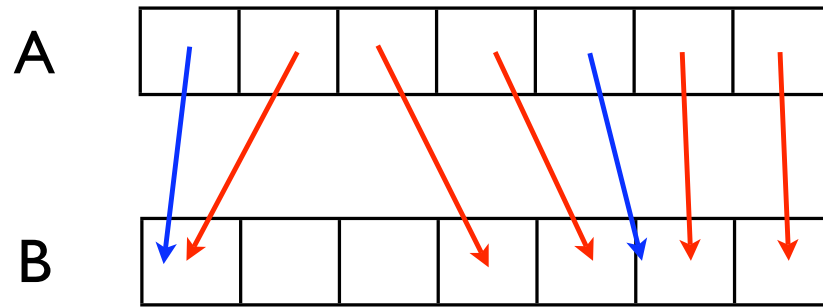

Warp pulls B into space of A, but is not dense in B.

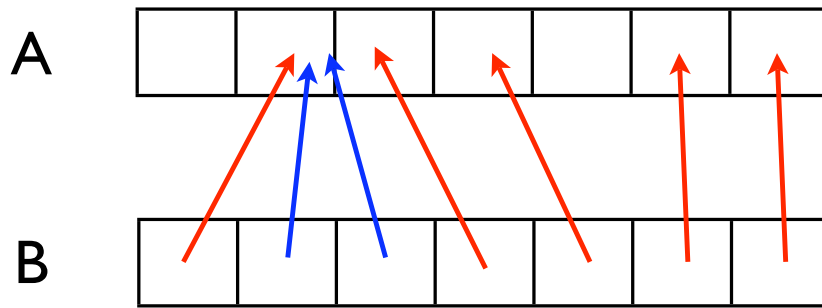


InverseWarp pulls A into space of B, but is not dense in A.

↑ Has inverse arrow  
in target domain

↑ No inverse arrow  
in target domain

Figure 6: Digital invertibility presents some limitations. Here, we see that invertibility is not exact but is gained only by interpolation. Thus, in three-dimensional scenarios in particular, there are fundamental limits to the degree of invertibility that may be achieved. The second and third voxels – from left – in image A undergo an expansion by a factor of 2. That is, under the mapping, 2 voxels are mapped to 4. This gives the definition of the Jacobian – computed by ANTSJacobian – which is a unitless measure defined by the ratio of volumes. Thus,  $J(\mathbf{x}) = V(\phi(\mathbf{x}))/V(\mathbf{x})$  where  $V$  represents the volume operation and  $\mathbf{x}$ , here, may be a small object. Thus, if  $\phi$  – the mapping – causes expansion, then  $J(\mathbf{x}) > 1$ .

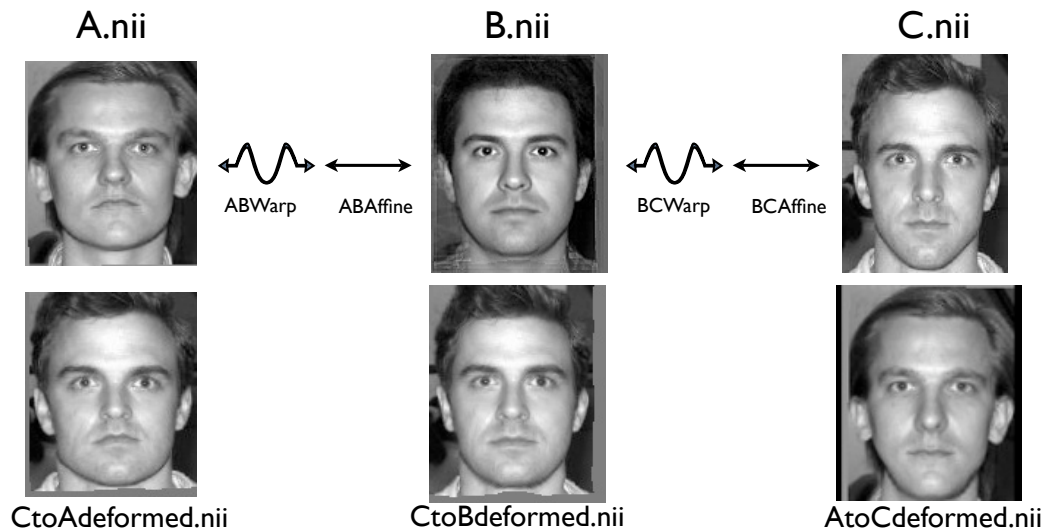

Read Warp Parameters left to right and pass to WarpImageMultiTransform - in short hand:

```
Warp C.nii CtoAdeformed.nii -R A.nii ABWarp ABAffine BCWarp BCAffine
```

Read right to left and take inverses to reverse the direction:

```
Warp A.nii AtoCdeformed.nii -R C.nii -i BCAffine BCInverseWarp -i ABAffine ABInverseWarp
```

Figure 7: Here, we detail how one applies a deformation and associated inverses. We also see how WarpImageMultiTransform can be used to concatenate a series of transformations.

**Diffeomorphic Transformations.** The elastic mapping space – as indicated above – may prove inadequate for some large deformation mapping scenarios. Figures 4 illustrate the changing performance one may get in switching from affine to elastic to diffeomorphic normalization. Figure 8 shows how one may use ANTS to achieve a state-of-the-art diffeomorphic mapping. The ANTS diffeomorphic model chosen for this example – symmetric normalization [10, 22] – is invariant to the order of the input images (although the affine mapping is not). An additional advantage of the diffeomorphic model over the elastic model is that both forward and inverse transformations are computed thus allowing one to deform fixed to moving and also moving to fixed. The transformation model chosen here – SyN[0.25] – may be replaced with other diffeomorphic transformation models. The step-size (0.25) impacts accuracy. Smaller is generally more accurate but takes more time to compute and may not capture as much deformation if the optimization is caught in a local minimum. The most general is global-in-time SyN via SyN[0.25, 2, 0.01], where the time step for integration is 0.01 (lower is more accurate). A fast approximate solution may be gained through exponential mapping via Exp[0.25, 10], where 10 integrations points are used. Update schemes for diffeomorphic optimization are similar to the elastic case and are described elsewhere [10, ?].


## 2.3 ANTS Similarity Terms and Other Functions

Here is a script that will let you experiment with different similarity term and transformation model combinations. A few are listed, with parameters that are a reasonable starting point.

```
# The user should define two images. here we use ANTS reference data
II=r16slice.nii.gz
```

## Symmetric Diffeomorphic Mapping with ANTS

Affine+ Diffeomorphic SyN

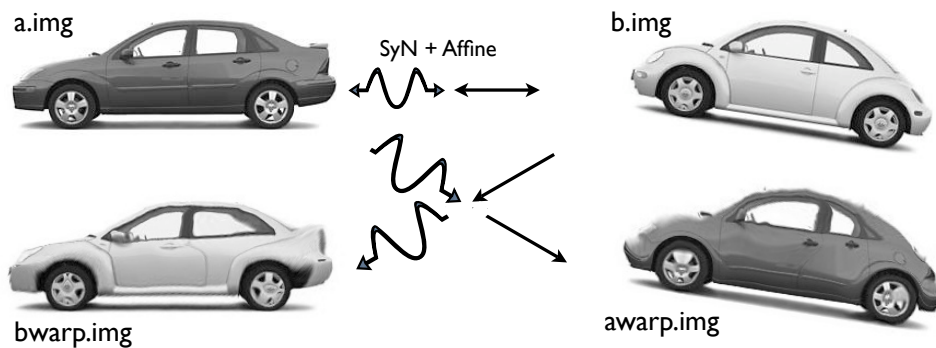

```
ANTS 2 -m PR[a.img,b.img,1,2] -i 100x100x10 -o ab.nii -t SyN[0.25] -r Gauss[3,0]
```

```
WarpImageMultiTransform 2 b.img bwarp.img -R a.img abWarp.nii abAffine.txt
WarpImageMultiTransform 2 a.img awarp.img -R b.img -i abAffine.txt abInverseWarp.nii
```

**Outputs:** abAffine.txt (invertible and composable with the warps)  
 abWarp.nii => abWarpxvec.nii & abWarpyvec.nii  
 abInverseWarp.nii => abInverseWarpxvec.nii & abInverseWarpyvec.nii

Figure 8: This example shows the benefit of the symmetric normalization model – invertibility, symmetry, highly deformable and accurate registration. The transformation model used here is highlighted in red.

```

JJ=r64slice.nii.gz
#
OUTPUTNAME=OUT # define the output prefix
ITS=" -i 100x100x30 " # 3 optimization levels

# different transformation models you can choose
TSYNWITHTIME=" -t SyN[0.25,5,0.01] -r Gauss[3,0] " # spatiotemporal (full) diffeomorphism
TGREEDYSYN=" -t SyN[0.15] -r Gauss[3,0] " # fast symmetric normalization
TELAST=" -t Elast[1] -r Gauss[0.5,3] " # elastic
TEXP=" -t Exp[0.5,10] -r Gauss[0.5,3] " # exponential

# different metric choices for the user
INTMSQ=" -m MSQ[ ${II}, ${JJ}, 1, 0] "
INTMI=" -m MI[ ${II}, ${JJ}, 1, 32] "
INTCC=" -m CC[ ${II}, ${JJ}, 1, 4] "
#
#
# these are the forward and backward warps.
#
IN VW=" -i ${OUTPUTNAME}Affine.txt ${OUTPUTNAME}InverseWarp.nii.gz "
FWDW=" ${OUTPUTNAME}Warp.nii.gz ${OUTPUTNAME}Affine.txt "
#
INT=$INTCC # choose a metric ( here , cross-correlation )
TRAN=$TGREEDYSYN # choose a transformation
if [[ $OPT == "Reg" ]] ; then

# run the registration
${ANTSPATH}/ANTS $DIM -o $OUTPUTNAME $ITS $TRAN $INT
# this is how you apply the output transformation
${ANTSPATH}/WarpImageMultiTransform $DIM ${II} ${OUTPUTNAME}IItoJJ.nii.gz -R ${JJ} $IN VW
${ANTSPATH}/WarpImageMultiTransform $DIM ${JJ} ${OUTPUTNAME}JJtoII.nii.gz -R ${II} $FWDW

elif [[ $OPT == "RegSeg" ]] ; then

# try something different -- segment the images and use them in ANTS
# to run multivariate registration.
IISEG=${OUTPUTNAME}A_seg.nii.gz
JJSEG=${OUTPUTNAME}B_seg.nii.gz
#
# get masks for atropos
#
${ANTSPATH}/ThresholdImage $DIM $II $IISEG 1 1.e9
${ANTSPATH}/ThresholdImage $DIM $JJ $JJSEG 1 1.e9
AtroposParams=" -d $DIM -m [0.1,1x1] -c [5,0] -i kmeans[3] "
${ANTSPATH}/Atropos $AtroposParams -a $II -o $IISEG -x $IISEG
${ANTSPATH}/Atropos $AtroposParams -a $JJ -o $JJSEG -x $JJSEG
# compute some segmentations and use them in labelguided mapping
LABELGUIDED=" -m PSE[ ${II}, ${JJ}, ${IISEG}, ${JJSEG}, 0.75, 0.1, 25, 0, 10] "

```

```
#
${ANTSPATH}/ANTS $DIM -o ${OUTPUTNAME} $ITS $TRAN $INT $LABELGUIDED
${ANTSPATH}/WarpImageMultiTransform $DIM ${II} ${OUTPUTNAME}IItoJJ.nii.gz -R ${JJ} $INVW
${ANTSPATH}/WarpImageMultiTransform $DIM ${JJ} ${OUTPUTNAME}JJtoII.nii.gz -R ${II} $FWDW

elif [[ $OPT == "BTP" ]] ; then
#
# now build a template from the two input images .
# in a real application, you would modify the call to btp
# such that it's more appropriate for your data.
# call buildtemplateparallel.sh -h to get usage.
#
${ANTSPATH}buildtemplateparallel.sh -d 2 -o ${OUTPUTNAME}BTP -c 0 $II $JJ
fi
```

Convergence occurs under two conditions: (1) either the maximum number of iterations are reached at a given optimization level or (2) the slope of the change in the optimization objective is negative or very small.

## 2.4 Choosing a Metric

ANTS supports both volumetric registration and point set registration. The image / point set similarity metrics in ANTS are unified in the form of a function on the images or the point sets:

**Similarity[fixedImage, movingImage, weight, parameters].**

The similarity type for the deformation transformation is specified by the `-m` option, which contains two parts: similarity type and the parameters inside the brackets. (Note: no white spaces exist between parameters.) The possible similarity metrics for volumetric images are:

- **Cross correlation estimate:** `"-m CC[fixedImage,movingImage,weight,radius]"`. This metric works well for intra-modality and some inter-modality image registration. For example, `-m CC[fixed.nii,moving.nii,1,5]` specifies:
  - the fixed image: `fixed.nii`
  - the moving image: `moving.nii`
  - weight for this metric is 1 (the weights are relative to the weights on the  $N$  other metrics passed to ANTS —  $N$  is unlimited). So, the weight,  $w_i$  on the  $i^{th}$  metric will be  $w_i = w_i / (\sum_i w_i)$ .
  - the region radius for computing cross correlation is 5
- **Mutual information:** `-m MI[fixedImage,movingImage,weight,number-of-histogram-bins]`. This metric works both well for intra-modality and inter-modality image registration. For example, the first three parameters in `-m MI[fixed.nii,moving.nii,1,32]` is similar to the example above in cross correlation, except that the last parameter means that the number of bins in computing mutual information is 32.
- **PR:** `-m PR[fixedImage,movingImage,weight,radius]`. This metric is deprecated. It still “works” but CC is preferred.



- **Mean square difference:** `-m MSQ[fixedImage,movingImage,weight,0]`. This metric works for intra-modality image registration. The last parameter 0 is a padding value of no real meaning. For example, `-m MSQ[fixed.nii,moving.nii,1,0]`.

ANTs also support registration of point sets. The supported formats for point sets can be found in I/O section. The similarity metrics for point sets are:

- **Point set expectation:**

```
-m PSE [fixedImage,movingImage,fixedPoints,movingPoints,weight,
pointSetPercentage,pointSetSigma,boundaryPointsOnly,
kNeighborhood,PartialMatchingIterations=100000]
```

- `fixedImage`: defines the space domain of the fixed point set.
- `movingImage`: defines the space domain of the moving point set.
- `fixedPoints/Image`: defines the coordinates of the fixed point set or label image. It can be an image with discrete positive labels, a VTK format point set file, or a text file. Details can be found in I/O section (TODO).
- `movingPoints/Image`: defines the coordinates of the moving point set or label image.
- `weight`: weight for this metric. 1 weights are relative to the weights on the  $N$  other metrics passed to ANTS —  $N$  is unlimited.
- `pointSetPercentage`: the percentage of points to be randomly sampled used in the registration.
- `pointSetSigma`: the standard deviation of the Parzen window used to estimate the expectation.
- `boundaryPointsOnly`: 1 (or “true”) means only the boundary points in the label image is used to drive registration.
- `kNeighborhood` is a positive discrete number. The first  $k$  neighbors are used to compute the deformation during the registration.
- `PartialMatchingIterations` controls the symmetry in the matching. This option assumes the complete labeling is in the first set of label parameters ... more iterations leads to more symmetry in the matching - 0 iterations means full asymmetry

- **Jensen-Tsallis BSpline**

```
-m JTB[fixedImage,movingImage,fixedPoints,movingPoints,weight,
pointSetPercentage,pointSetSigma,boundaryPointsOnly,kNeighborhood,alpha,
meshResolution,splineOrder,numberOfLevels,useAnisotropicCovariances]
```

- `fixedImage`: defines the space domain of the fixed point set.
- `movingImage`: defines the space domain of the moving point set.
- `fixedPoints`: defines the coordinates of the fixed point set. It can be an image with discrete positive labels, a VTK format point set file, or a text file. Details can be found in I/O section (TODO).
- `movingPoints`: defines the coordinates of the moving point set.

fixed image | moving image | MSE

Figure 9: registration using mean square intensity difference

- weight: weight for this metric. The weights are relative to the weights on the  $N$  other metrics passed to ANTS —  $N$  is unlimited.
- pointSetPercentage: the percentage of points to be randomly sampled used in the registration.
- pointSetSigma: the sigma for the Parzen window used to estimate probabilities.
- boundaryPointsOnly: [TODO] 1 (or “true”) means only the boundary points in the point sets are used to drive registration.
- kNeighborhood is a positive discrete number. The first  $k$  neighbors are used to compute the deformation during the registration.
- alpha
- meshResolution
- splineOrder
- numberOfLevels
- useAnisotropicCovariances

ANTS has three affine registration similarity metrics for volumetric images, which are specified using `--affine-metric-type`:

- Mutual information, specified as `--affine-metric-type MI`. This usually works for both inter and intra-modalities in 3D. Also, the options used in computing mutual information in affine registration can be specified as `--MI-option N1xN2`. The first parameter  $N1$  specifies the number of bins. The second parameter  $N2$  specifies the number of samples. For example: `--MI-option 32x8000`.
- Cross-correlation `--affine-metric-type CC`. CC has a wide capture range.
- Mean square difference, specified using `--affine-metric-type MSE`. There is no options for this metric.

Fig. 2.4 shows the registration result using intensity difference with the following command:

```
ANTS 2 -m MSQ[r16slice.nii,r64slice.nii,1,0] -r Gauss[3,0] -t SyN[0.5] -i 50x50x30
```

Here is also an example script to register a pair of images using mean square intensity difference and computing the metrics of the registration image.

```
#use intensity difference with radius 0 -- radius no effect on intensity difference
ANTS 2 -m MSQ[r16slice.nii,r64slice.nii,1,0] -r Gauss[3,0] -t SyN[0.5] -i 50x50x30
WarpImageMultiTransform 2 r64slice.nii resMSQ.nii Warp.nii Affine.txt
MeasureImageSimilarity 2 0 r16slice.nii resMSQ.nii metricexamplelog.txt
MeasureImageSimilarity 2 1 r16slice.nii resMSQ.nii metricexamplelog.txt
MeasureImageSimilarity 2 2 r16slice.nii resMSQ.nii metricexamplelog.txt
ConvertToJpg resMSQ.nii resMSQ.jpg
```

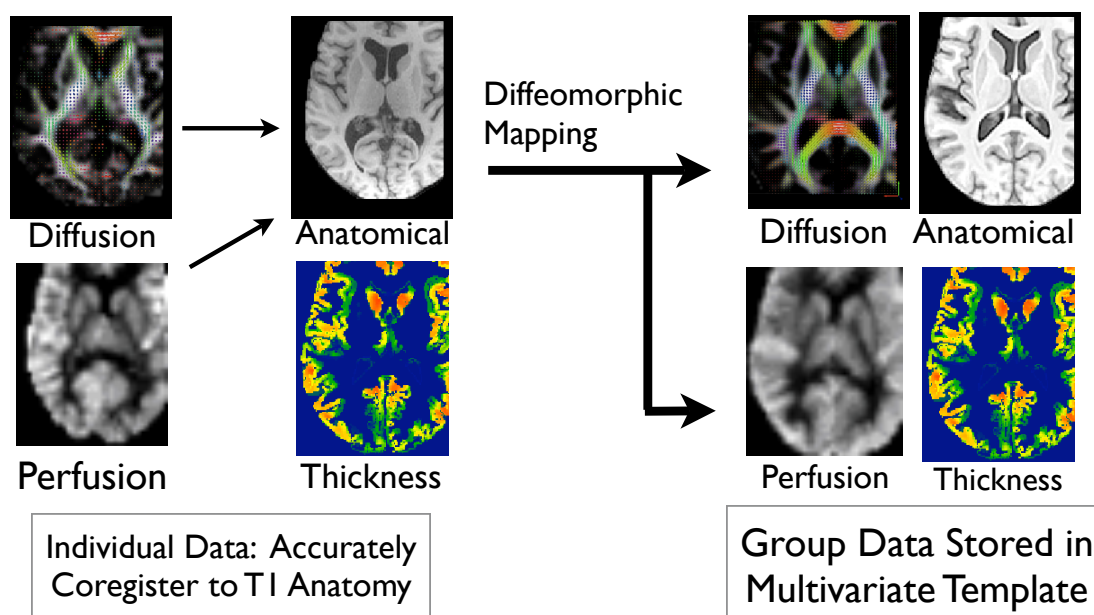


Figure 11: The current version of ANTS may perform normalization of multiple modalities by combining *intrasubject*, *intermodality* mappings with *intersubject*, *intramodality* maps to a group template. In brain imaging, the intersubject maps are usually guided by the T1 component.

## 2.5 Diffusion Tensor Normalization

The ImageMath program – via TensorFA and TensorMeanDiffusion – may derive scalar images from DTI. Such images may be used to distortion correct DTI to T1 or to map DTI together. See the ANTS/Scripts program called `antswithdt.sh` for a ready to go example. Figure 10 shows what might happen if your tensor entries are stored in the wrong order. ANTS expects the nifti standard ordering of the DT six vector.

## 2.6 Multivariate Normalization with ANTS

Multivariate normalization may be performed in two ways with ANTS. First, as shown in figure 11, intra-subject mappings may be used to conform all modalities into a common subject space defined by the “anchor” modality. In brain imaging, this is usually the T1 component.

A second type of multivariate normalization is able to use the T1 and DTI components directly in the optimization of the mapping. An example of this type of call is below, where we assume that one has resampled the DT data to the space of T1:

```
ANTS 3 -m CC[T1template.nii,T1subject.nii,1.25,4]
      -m CC[FAtemplate.nii,FAsubject.nii,1,4]
-o MultiVar -i 30x30x20 -r Gauss[3,0] -t SyN[0.25]
```

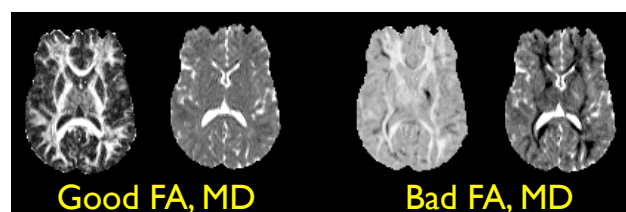


Figure 10: The FA and mean diffusion of the tensor may be corrupted, as at right, if the order of the DT entries is wrong.

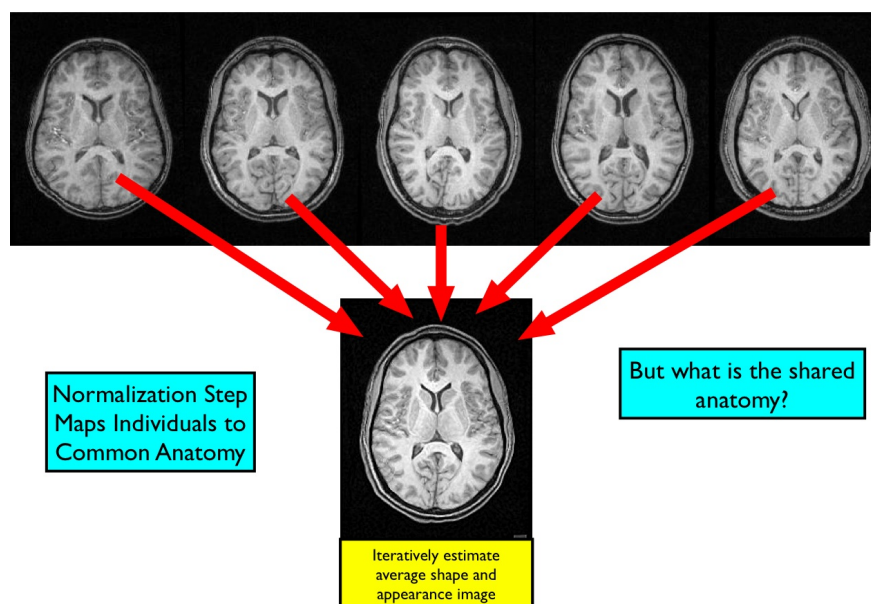


Figure 12: This example may be recreated by the user. The shared anatomy – across this dataset – is recovered in the derived optimal template.

In this case, dense features within white matter are gained by using the DT component in addition to the T1 component. The T1 CC metric, here, is weighted slightly more than the FA CC metric, 1.25 vs. 1.0. *The convergence criterion for the multivariate scenario is a slave to the **last** metric you pass on the ANTS command line.*

## 2.7 Optimal Template Construction with ANTS

A useful script for setting up optimal template construction is available in ANTS/Scripts. This method is described in [6, 16, 22, 30, 9], where more recent publications are more descriptive and incorporate more recent features and optimization strategies. Two versions of this algorithm are available : serial and parallel. Data and a script for building an average face template is available the ANTS/Examples/Data directory. Two key points about optimal templates: 1. The outcome stabilizes at around ten images, for most populations. That is, if we randomly select images from individuals within the same demographic distribution, we will end up with very similar average images. 2. Optimality, for the ANTS-SyN approach, is defined by the minimum shape and appearance distance image. See the template page for some examples of this and users may contact ANTS developers for previously derived template images that may be useful.

### A Concrete Template Construction Example, Step-by-Step

For this example, we will use the data in "ANTS/Examples/Data/B\*.tiff." The files are B1.tiff, B2.tiff, B3.tiff, B4.tiff and B5.tiff. They are single slices from a dataset of real MRI brain images – individuals around 15 years of age. Our goal is to derive a "most representative" single image from this population. The steps are:

1. Make a directory for this example. Copy or link all the B\*.tiff images into it.

2. Copy the file `ANTS/Scripts/buildtemplateparallel.sh` into that directory.
3. On the command line, within that directory, run the following command: `bash buildtemplateparallel.sh -h` to get usage.
4. Following the usage, you can write a script that contains a call like this:  
`buildtemplateparallel.sh -d 2 -o BTP -c 0 B*tiff`
5. Call the script in step 4 and then wait for the output. The template will be called (in this example) `BTPtemplate.nii.gz`.
6. If the script does not produce the desired output, check that you defined `$ANTSPATH` in your path. You may also check and or define the path to the ANTs programs explicitly by editing this variable in `buildtemplateparallel` and `antsIntroduction`.

The script offers an option ( `-c 1` for `qsub` , `-c 2` to attempt to use multiple cores on your machine ) so that you can send the whole thing to distributed computing. Typically, one would use `voxb` or the `qsub` program available in the Sun Grid Engine. Figure 12 shows the data and a result that may be recreated by the user.

## 2.8 Notes on Large Deformation Mapping

Figure 13 defines what one might expect from a high-resolution, large deformation, successful normalization. Major and many minor features are well aligned between these images to an extent that is approaching the limits allowable by biological variation. **Turning "Failures" into Successes:** Below are some pointers to follow if you are unable to recreate such normalization quality. Usually, reasons for registration failure are one of a very few things:

1. The initialization is so off mark that affine registration fails – thus meaning all subsequent deformable registration will not be meaningful.
2. The information within headers is inconsistent e.g. origins are defined in different ways, the "directions" are not correct or some combination of these. The `PrintHeader` executable can help one in debugging this type of problem. Also, the `ImageMath` function `CompareHeadersAndImages` can fix some problems of this type.
3. The similarity or transformation model is inappropriate or has too small a capture region for the problem.

An example of this final point is given here. A recent study used the following call to ANTS across a dataset of elderly subjects:

```
ANTS 3 -m CC[template.nii.gz,subject.nii.gz,1,2] -i 10x50x50x20
-o subjectmap.nii.gz -t SyN[0.25] -r Gauss[3,0]
```

This succeeded for all subjects but one. This subject had severe neurodegeneration which caused "under-normalization" to result. A larger deformation mapping was required and so we increased the maximum number of iterations (the "-i" parameter). One might also increase the gradient-descent step-size (`SyN[0.25]` → `SyN[0.5]`). Large deformation mapping is challenging not only because of the amount of deformation, but also because the "true" solution is difficult to find. Keeping this in mind, we also increased the "span" (radius) of the correlation window in the similarity metric (which also increases computation time). These modified parameters succeeded:

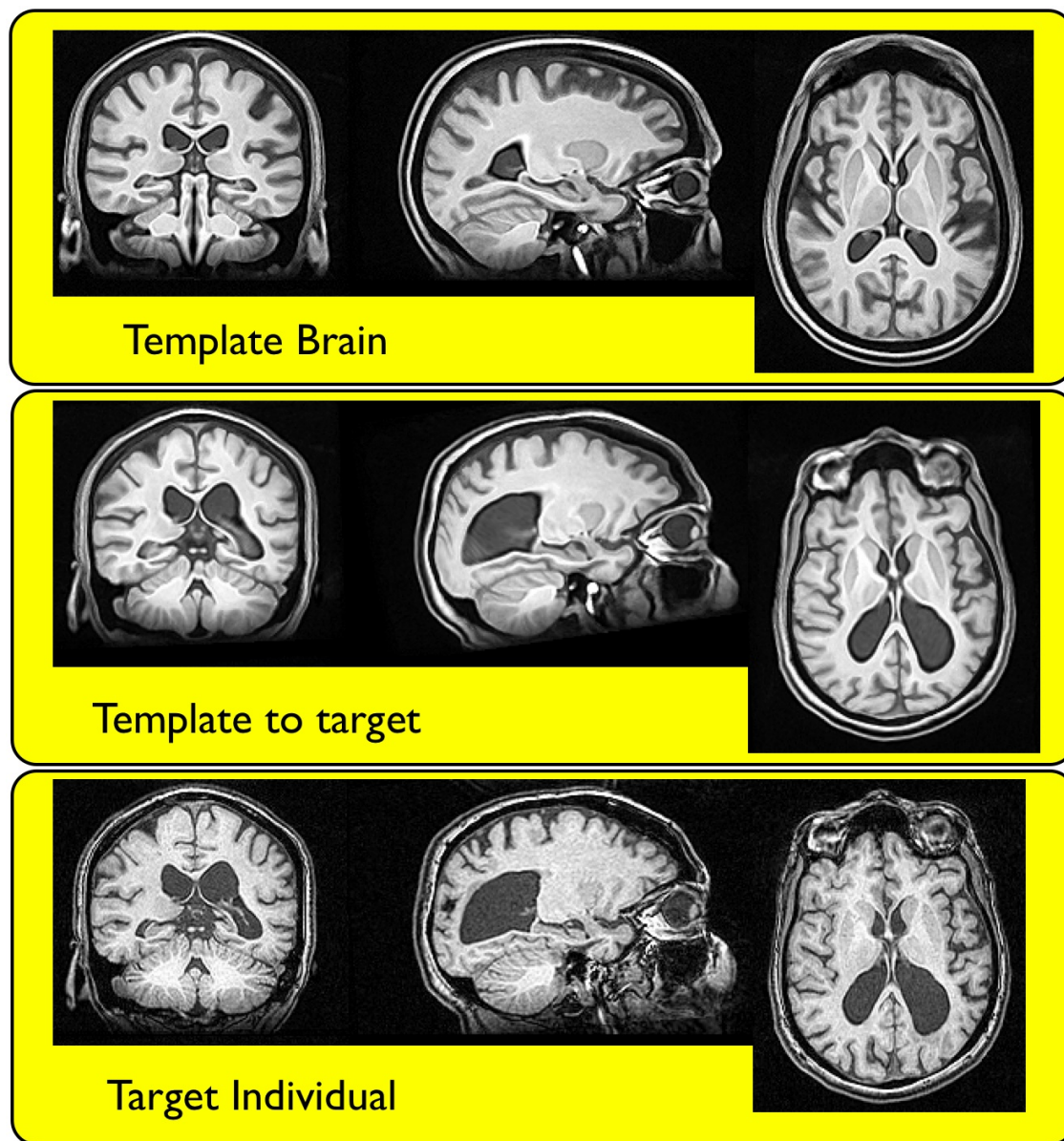


Figure 13: ANTS succeeds in a challenging normalization scenario. Comparing the template mapping to the individual shows that much of the cortex is well-aligned, as is the hippocampus, despite the relatively large difference between the initial template and the target. Here, one might note a limitation of whole brain mapping: occipital lobe sulcal variation is highly idiosyncratic and extremely difficult to capture, in particular when there is such severe atrophy nearby.



```
ANTS 3 -m CC[template.nii.gz,subject.nii.gz,1,4] -i 100x100x100x20  
-o subjectmap.nii.gz -t SyN[0.25] -r Gauss[3,0]
```

Key changes were to increase the radius of the correlation and to allow more iterations at the coarsest resolution. This two coarsets levels of the computation take only about 10 percent of the time (about 5-10 minutes depending on the machine) but accounts for the large majority of the shape variation. Note that this is a 3D example – that is why the eyes appear only in the frame at left.

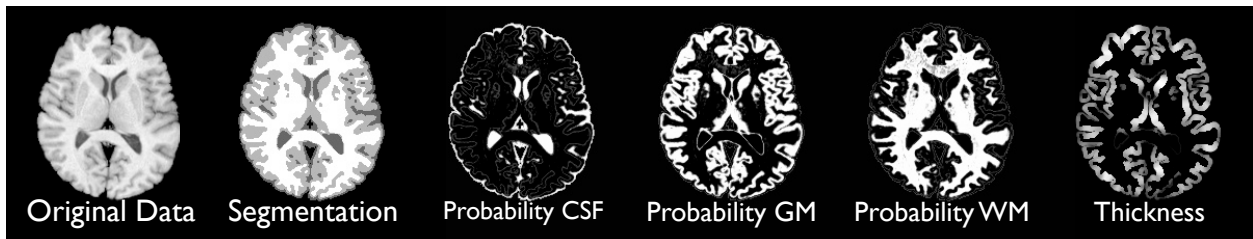


Figure 14: ANTS provides basic Markov Random Field regularized Gaussian model tissue segmentation. The ANTS program `N3BiasFieldCorrection` is very valuable as preprocessing for this naive approach to segmentation. The last panel – far right – shows the thickness derived from the white matter and gray matter probabilities where a prior on thickness was used to prevent thickness overestimation. Thickness was derived with DiReCT [17], a software tool available in binary form only – a future release may provide code.

### 3 Image Segmentation and Labeling

ANTS has tools for both tissue based segmentation and prior-based segmentation that leverages spatial priors, usually based on a template mapping.

#### 3.1 Basic Segmentation

An advanced MRF segmentation algorithm is available in ANTS via the `Atropos` tool, which is still under development. We apply the segmentation to the example data in: `ANTS/Examples/Data/r16slice.nii`.

```
ThresholdImage 2 r16slice.nii mask.nii.gz 0.001 Inf
Atropos -d 2 -a r16slice.nii -o [r16slice_seg.nii.gz,r16slice_prob_%02d.nii.gz] \
  -c [10,0.001] -m [0.1,1x1x1] -i kmeans[3] -x mask.nii.gz
```

The parameter `kmeans[3]` indicates 3 tissues plus background are sought. The `[0.05,1]` indicates that some MRF smoothing will be applied to the labels. Other options are defined in the `Atropos` help. Output is shown in figure 14 and – with priors – in figure 15. The document `atropos.pdf` in the Documentation directory provides much more explanation.

#### 3.2 Prior and Template-based Image Segmentation

The same algorithm may be augmented to perform Prior-based segmentation.

```
Atropos -d 2 -a r16slice.nii -o [r16slice_seg.nii.gz,r16slice_prob_%02d.nii.gz] \
  -c [10,0.001] -m [0.05,1] \
  -i PriorProbabilityImages[3,r16slice_prob_%02d.nii.gz,0.5] -x mask.nii.gz
```

The `PriorProbabilityImages` parameter 3 indicates 3 tissues plus background are sought. The `PriorProbabilityImages` term 0.5 indicates we weight the priors equally as the data term and use a spatially varying set of Gaussians to estimate the segmentation. In this way, data from the priors may be used to modify and guide the segmentation in a locally varying way, accommodating for both inhomogeneity and the different imaging signature that different tissues provide. Prior images should be of the same size and dimension as the input data and should have intensities in the range  $[0, 1]$ .



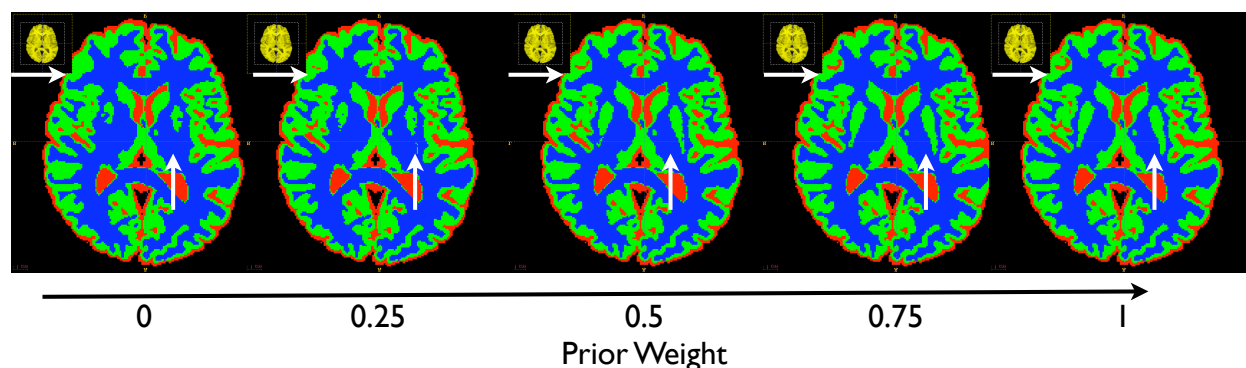


Figure 15: The ANTs program *Atropos* extends basic Markov Random Field regularized Gaussian model tissue segmentation to include priors, which allow spatially varying tissue models to determine segmentation. The ANTs program *N4BiasFieldCorrection* is critical as preprocessing when using this approach to segmentation. Here, we see the improvement in segmentation as the locally varying prior models are weighted more heavily. The last panel – far right – shows the segmentation derived from data driven with an initialization founded purely on prior models based on template mapping. The prior model brings out the caudate and some CSF – highlighted by arrows – as their weight increases.

### 3.3 Cortical Thickness

Two forms of cortical thickness estimation from probability maps are available in ANTS: first, the traditional Laplacian cortical thickness estimation and, second, the more recently developed Diffeomorphic Registration-based Cortical Thickness (DiReCT) [17]. Both methods estimate thickness of an object from probability or binary images that define object boundaries. This tool is mainly of interest in brain mapping and cardiac imaging for morphometric studies. Cortical thickness, for instance, is known to correlate with language development and IQ in adolescents. The laplacian approach to cortical thickness may be used via the program *LaplacianThickness*, which implements the method described in “Three-dimensional mapping of cortical thickness using Laplace’s Equation” by Jones, et al [1]. Super-sampling and controlling segmentation accuracy for input to this program is up to the user. Otherwise, the Laplacian method may grossly overestimate thickness in closed sulci.

### 3.4 User/Label-Guided Normalization for Hippocampus Mapping

See the link for details on hippocampus labeling: <http://picsl.upenn.edu/ANTS/hipptutorial.php>. **Expectation-Based Point-Set Registration.** Here, we apply the expectation-based point set registration method for mapping labeled points sets, as described in [25]. ITK-SNAP may be used to label images and exported segmentation images may be input to the PSE metric below, as labeled data. The Frown and Smile data is used as example. This data is available in the ANTS/Examples/Data/ directory.

```
ANTS 2 -o PSEtest -i 91x70x55x40x30 -r Gauss[3,0.] -t SyN[0.2]
      -m PSE[Frown.nii,Smile.nii,Frown.nii,Smile.nii,0.75,0.1,11,0,10]
      -m MSQ[Frown.nii,Smile.nii,1,0] --number-of-affine-iterations 0
```

```
WarpImageMultiTransform 2 Frown.nii FrownToSmile.nii -R Smile.nii
-i PSEtestAffine.txt PSEtestInverseWarp.nii
```

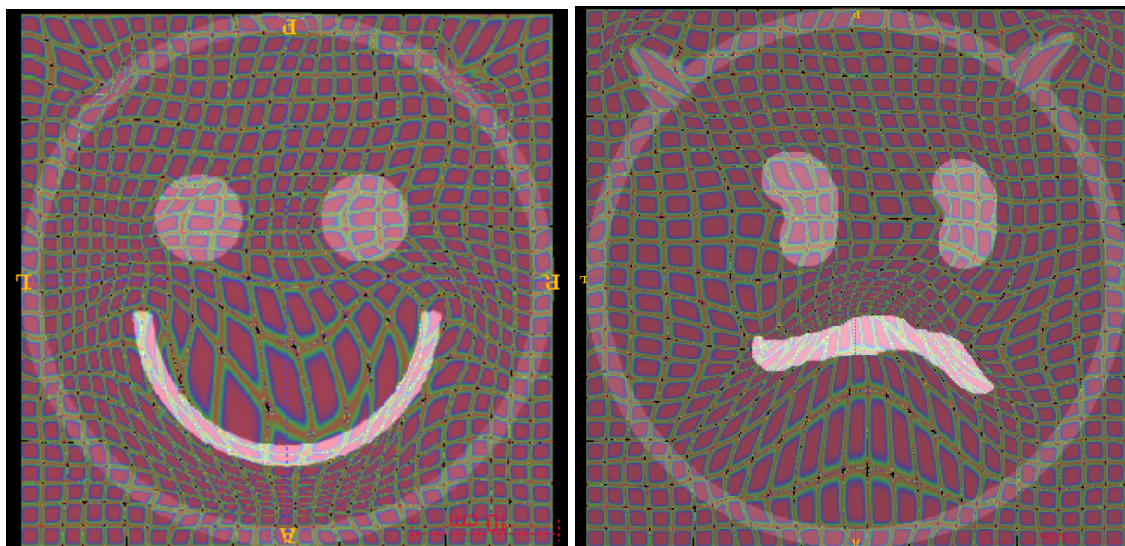


Figure 16: Expected output for the frown to smile shows a smooth, though large deformation. The grids are overlaid on the deformed images.

```
WarpImageMultiTransform 2 Smile.nii SmileToFrown.nii -R Frown.nii  
PSEtestWarp.nii PSEtestAffine.txt
```

```
CreateWarpedGridImage 2 PSEtestInverseWarp.nii grid1.nii  
CreateWarpedGridImage 2 PSEtestWarp.nii grid2.nii
```

This example should run on the downloaded ANTS data so you may see the results.

## 4 Application to Studies

### 4.1 Brain Mapping in the Presence of Lesions

Many cases violate the basic assumption of a one-to-one mapping from a template image to a target image. Brain lesions caused by stroke or traumatic brain injury are a common instance of this.

The ANTS toolkit handles this situation by a constrained cost-function masking approach. In short, the mapping within an excluded region (as determined by an inclusive mask) is determined by the solution outside of the region. To achieve this with ANTS, one would use a call similar to this:

```
ANTS 3 -m CC[tp22_s1.nii,template.nii.gz,1,4] -i 50x20x0
      -o tp22map -t SyN[0.25] -x mask.nii.gz -r Gauss[3,0]
```

The additional parameter needed for this approach is the `-x mask.nii.gz` which is the mask defined in the "fixed" (or individual) image space. Mask values above 0.5 are considered non-outlier or non-lesioned regions. The labeling may be performed with ITK-SNAP. Figure 17 illustrates the approach. Once the ANTS mapping is solved, one is able to estimate the location of the lesion in the template space by mapping the inverse lesion mask to the template. To take the negative image of the inclusive mask and warp the mask to the template :

```
ImageMath 3 lesionmask.nii.gz Neg inclusivemask.nii.gz
WarpImageMultiTransform 3 lesionmask.nii.gz lesiontotemplate.nii.gz
      -R template.nii.gz -i tp22mapAffine.txt tp22mapInverseWarp.nii
```

More details on this approach are being prepared for publication. See also: [http://www.ncrrn.org/papers/methodology\\_papers/sp\\_norm\\_kim.pdf](http://www.ncrrn.org/papers/methodology_papers/sp_norm_kim.pdf).

### 4.2 Statistical Mapping with ANTS: Morphometry, Function, Jacobian, Thickness

ANTS has been applied in a wide array of imaging studies [5, 4, 14, 17, 23, 24, 25, 29, 30, 12, 10, 20, 22, 27, 28, 31, 1, 3, 15]. All of these studies benefit in some way from normalization whether the topic is cortical thickness, Jacobian-based morphometry, volumetric morphometry or functional studies. From a broad perspective, each of these applications requires the same steps:

1. Preprocess images – bias correction (e.g. via the `N4BiasFieldCorrection` program), segmentation (via `Atropos`), etc. Potentially construct an optimal template.
2. Normalize images – run `ants.sh` to map a population of images to a template and store the deformations.
3. Derive any data from the deformation that may be necessary, e.g. the log-Jacobian, via `ANTSJacobian`.
4. Apply the warp to any images one may want to analyze in template space.
5. Perform statistics in template space over the region of interest, e.g. all cortex.

We now illustrate this procedure using ANTS tools and example images from `ANTS/Examples/Data`. The expected output is in figure 18.

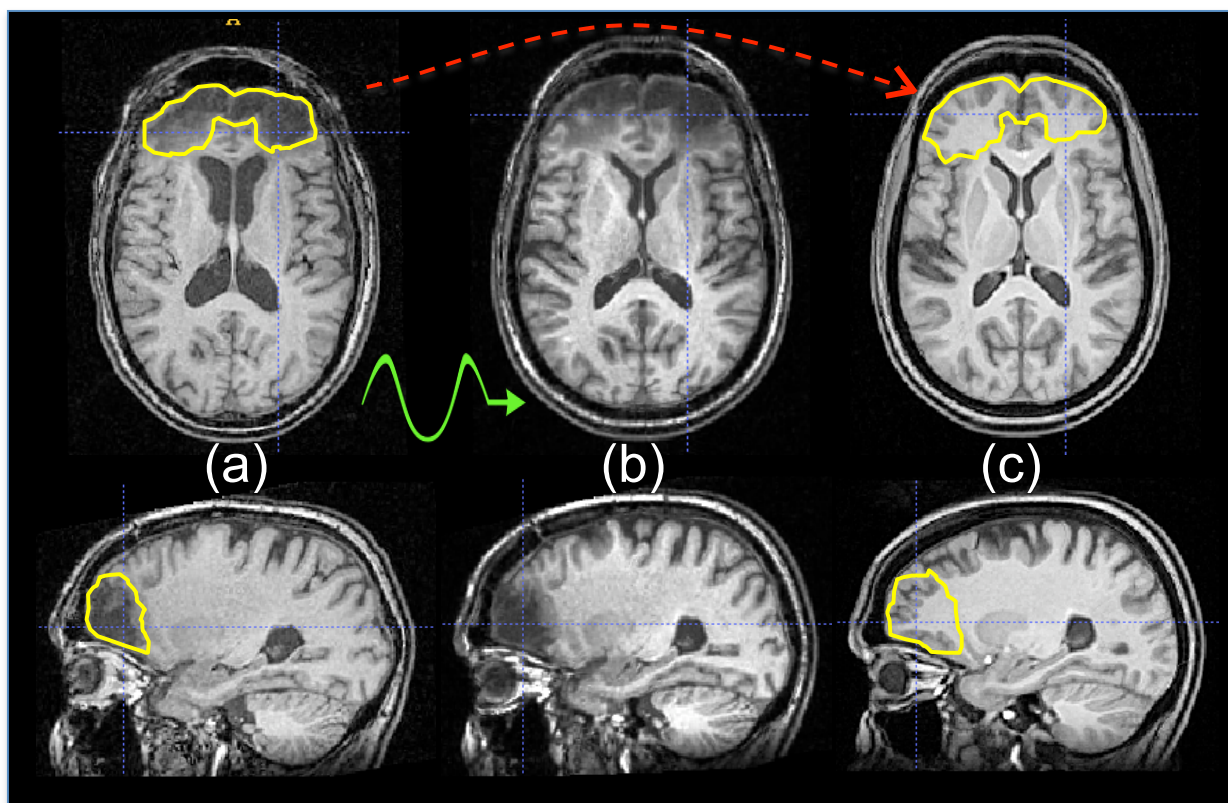


Figure 17: Here, we see a semi-automated approach for lesion analysis based on diffeomorphic normalization. The user outlines the lesion site in the subject space, as shown in (a). Diffeomorphic mapping is then used to deform the subject image in (a) into the space of the reference template. The deformed subject is shown in panel (b). Panel (c) shows the reference template space with the estimated position of the subject's lesion overlaid on the healthy tissue of the template. This approach enables one to compare the subject image against prior information stored in the template space such as anatomical labels or statistics on the appearance of lesioned and/or normal tissue.

```

# 0. Here, we use r16 as a template and r64 as a subject T1 image.
# 1. Segment the input images
ThresholdImage 2 r16slice.nii mask.nii.gz 0.001 Inf
Atropos -d 2 -a r16slice.nii -o [r16slice_seg.nii.gz,r16slice_prob_%02d.nii.gz] \
    -c [10,0.001] -m [0.05,1] -i kmeans[3] -x mask.nii.gz
ThresholdImage 2 r64slice.nii mask.nii.gz 0.001 Inf
Atropos -d 2 -a r64slice.nii -o [r64slice_seg.nii.gz,r64slice_prob_%02d.nii.gz] \
    -c [10,0.001] -m [0.05,1] -i kmeans[3] -x mask.nii.gz
# 2. Compute the deformation
sh ants.sh 2 r16slice_prob_02.nii.gz r64slice_prob_02.nii.gz SPM 100x100x100
# 3. Smooth the segmented image
SmoothImage 2 r64slice_prob_02.nii.gz 1. sr64slice_prob_02.nii.gz
# 3. Warp the segmented image by the
WarpImageMultiTransform 2 sr64slice_prob_02.nii.gz sr64slice_prob_02_norm.nii.gz \
    -R r16slice.nii SPMWarp.nii.gz SPMAffine.txt
# 4. Compute the log Jacobian from the transformation.
ANTSJacobian 2 SPMWarp.nii.gz SPM 1 # take the log
gzip SPMlogjacobian.nii
# 5. Multiply the log-jacobian by the gray matter probability.
ImageMath 2 SPMlogjacobianmask.nii.gz m SPMlogjacobian.nii.gz sr64slice_prob_02_norm.nii.gz
# 6. Estimate cortical thickness --- use svn revision 684 or greater (not in official release)
KellySlater 2 r64slice_seg.nii.gz r64slice_prob_03.nii.gz r64slice_prob_02.nii.gz \
    r64slice_thickness.nii.gz 0.02 50 8 0 0.5
# 7. Warp to template space.
WarpImageMultiTransform 2 r64slice_thickness.nii.gz r64slice_thickness_norm.nii.gz \
    -R r16slice.nii SPMWarp.nii.gz SPMAffine.txt
# 8. Repeat steps 1-7 for a population and run statistics.

```

This example uses a standard Jacobian-based approach to morphometry. Here, we use the log-Jacobian because it is symmetric about zero and mask with the gray matter segmentation to restrict the analysis to the cortex. The Jacobian is discussed in figure 6 and visualization of the Jacobian is shown in figure 19. This Jacobian comes from the mapping shown in figure 13. Similar analyses may be performed on thickness images, functional images or fractional anisotropy images derived from the diffusion tensor modality. In all these cases, steps similar to those above would be performed.

### 4.3 Statistics with ANTS and **R**

**R** is an open-source, cutting-edge statistical package used by statisticians world-wide. ANTS is designed to interface with **R** via an ImageMath tool that converts image sets into an **R** compatible format. Three steps are required. First, convert image data into vector and matrix formats through a mask of your statistical region of interest. Second, read that data into **R** and run your statistics. Third, write out the data from **R** and convert back to an image. Here are the steps, in code examples,

1. Create a vector for your output:

```
ImageMath 3 testresultvec.nii ConvertImageSetToMatrix 1 mask.nii image1.nii
```



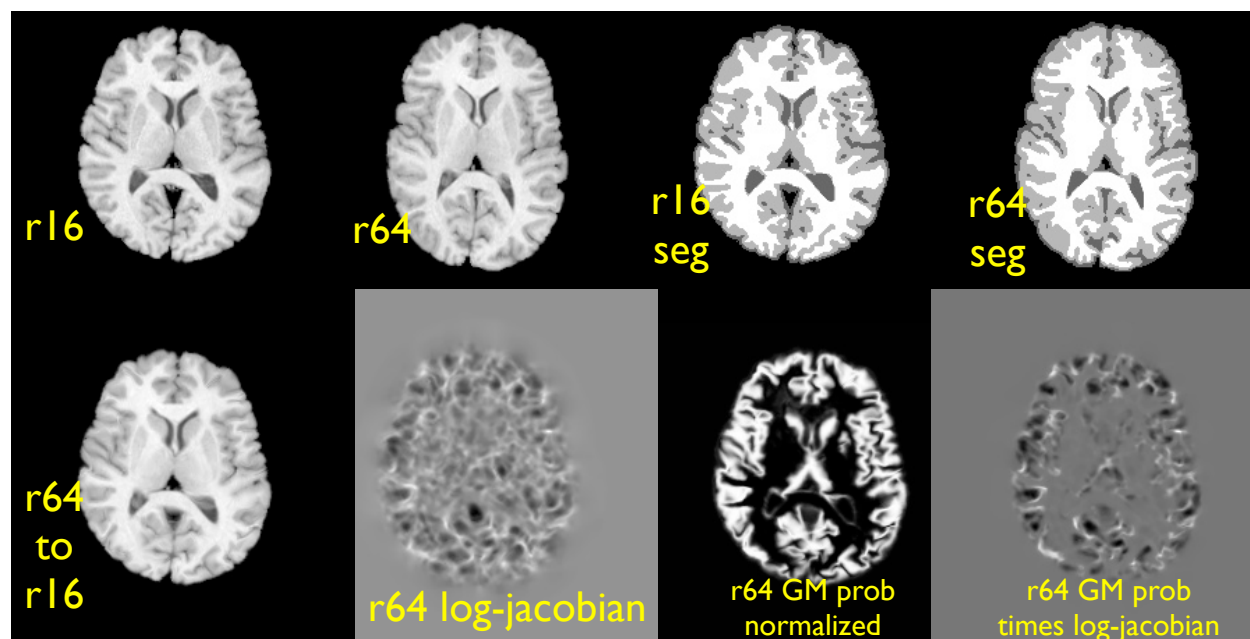


Figure 18: The output from the ANTS morphometry example.

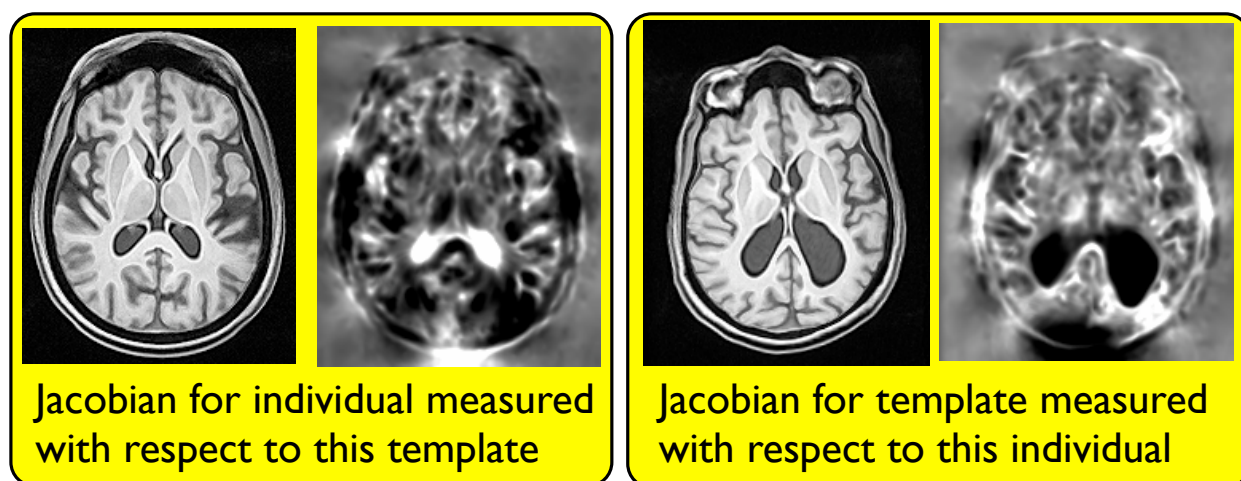


Figure 19: The Jacobian: Note that the bright values in the template image ventricles (left) indicate that the ventricles are relatively larger in the subject image. Similarly, the dark jacobian values in the individual image show that the template has smaller ventricles. This example is described on the large deformation page Large Def. Once one gains the Jacobian (or more appropriately, the log-Jacobian), then one may compute statistics across the population.

Here, we assume all images are called image1.nii, image2.nii, ... , imageN.nii so we can use a wildcard to pass to ImageMath and create a matrix for all your input data:

```
ImageMath 3 datamatrix.nii ConvertImageSetToMatrix 1 mask.nii image*.nii
```

2. Run **R** by reading your data above, performing whatever statistics you want and then writing out.
3. Convert the test result back to an image.

```
ImageMath 3 testresult.nii ConvertVectorToImage mask.nii testresultvec.nii
```

Most statistical requirements may be met with this setup.

## 5 Dependencies and Related Software

### 5.1 Dependencies and Compilation

ANTs depends on the most recent version of the Insight ToolKit (ITK) – [www.itk.org](http://www.itk.org). File types that may be read and written are all of those available in itk and a few more <http://picsl.upenn.edu/ANTs/ioants.php>.

ANTs and ITK are both compiled through pointing CMake – [www.cmake.org](http://www.cmake.org) – to a CMakeLists.txt file. The ANTs CMakeLists.txt is in the Examples directory.

See below for more details on compilation and binaries:

**Compile and download:** <http://picsl.upenn.edu/ANTs/download.php>. We maintain a dashboard to give users an idea of what to expect from compilation:

<http://picsl.upenn.edu/cdash/index.php?project=ANTs>

### 5.2 Visualization and Quantification of ANTs Results

Use **itk SNAP** : [www.itksnap.org](http://www.itksnap.org) or mricro / mricron for viewing images. The ANTs program StackSlices can be used to conveniently scan through a dataset in any of its axes. It is useful for checking both input data quality and output normalization quality. MeasureImageSimilarity will supply numerical measures of image similarity via three different metrics.

### 5.3 Statistics Beyond ANTs

Use SPM, **R** or nfm (from Chris Rorden's mricron) for computing population statistics.

### 5.4 Pipelining with ANTs

Use PipeDream for automating large-scale processing:

**PipeDream Homepage** : <https://sourceforge.net/projects/neuropipedream/>

PipeDream is integrated with ANTs and automates more complex studies such as multivariate diffusion tensor and T1 cortical thickness studies, longitudinal mapping and reconstruction of nifti images from Dicom.



## 6 Annotated Bibliography

The *de facto* ANTs reference is now [“A reproducible evaluation of ANTs similarity metric performance in brain image registration.”](#) (this is a link). The original statement of the symmetric normalization and template construction methodology was given in [6] with a more recent template paper here [9]. A follow up study that used landmark guidance to compare the chimpanzee cortex to the human cortex was published here [8] – this study used *in vivo* MRI and template-based normalization to confirm volumetric numbers derived from an early 20th century post-mortem study comparing one human and one chimp. This conference article has some additional detail and alternative updates to the methodology, in particular application to shape-based interpolation [7]. Network based studies were performed here [18, 19]. The main SyN paper is here [13]. Applications to neurodegeneration are here [2, 11, 21, 14, 17, 30, 24]. Hippocampus focused work is here [26, 30]. The main evaluation papers include [13] and [23] for the cortex and deep brain structures whereas [26] evaluates the use of automated and semi-automated normalization for high-throughput hippocampus morphometry. An additional evaluation paper is being developed.

## References

- [1] Geoffrey K Aguirre, Andres M Komeromy, Artur V Cideciyan, David H Brainard, Tomas S Aleman, Alejandro J Roman, Brian B Avants, James C Gee, Marc Korczykowski, William W Hauswirth, Gregory M Acland, Gustavo D Aguirre, and Samuel G Jacobson. Canine and human visual cortex intact and responsive despite early retinal blindness from rpe65 mutation. *PLoS Med*, 4(6):e230, Jun 2007. [4.2](#)
- [2] B. Avants. *Shape optimizing diffeomorphisms for medical image analysis*. PhD thesis, University of Pennsylvania, 2005. [6](#)
- [3] B. Avants, C. Anderson, and M. Grossman. Tauopathic longitudinal gray matter atrophy predicts declining verbal fluency: A symmetric normalization study. In *Human Brain Mapping*, 2007. [4.2](#)
- [4] B. Avants, P. A. Cook, J. T. Pluta, J. nad Duda, H. Rao, J. Giannetta, H. Hurt, S. Das, and J. Gee. Follow-up on long term effects of prenatal cocaine exposure/poly-substance abuse on the young adult brain. In *Pediatric Academic Society, Annual Meeting*, 2009. [4.2](#)
- [5] B. Avants, P. A. Cook, J. T. Pluta, J. nad Duda, H. Rao, J. Giannetta, H. Hurt, S. Das, and J. Gee. Multivariate diffeomorphic analysis of longitudinal increase in white matter directionality and decrease in cortical thickness between ages 14 and 18. In *Human Brain Mapping, 15th Annual Meeting, oral presentation.*, 2009. [4.2](#)
- [6] B. Avants, C. L. Epstein, and J. C. Gee. A method for conformally mapping simply connected domains to the unit disc. in preparation, 2004. [2.7](#), [6](#)
- [7] B. Avants, C. L. Epstein, and J. C. Gee. Geodesic image interpolation: Parameterizing and interpolating spatiotemporal images. In *ICCV Workshop on Variational and Level Set Methods*, pages 247–258, 2005. [6](#)
- [8] B. Avants, C. L. Epstein, and J. C. Gee. Geodesic image normalization in the space of diffeomorphisms. *Mathematical Foundations of Computational Anatomy*, pages 125–133, 2006. [6](#)
- [9] B. Avants, P. Yushkevich, J. Pluta, and J. C. Gee. The optimal template effect in studies of hippocampus in diseased populations. *Neuroimage*, page in press, 2009. [2.7](#), [6](#)

- [10] B. B. Avants, C. L. Epstein, M. Grossman, and J. C. Gee. Symmetric diffeomorphic image registration with cross-correlation: Evaluating automated labeling of elderly and neurodegenerative brain. *Med Image Anal*, 12(1):26–41, Feb 2008. [1](#), [2](#), [2.2](#), [4.2](#)
- [11] B. B. Avants, C. L. Epstein, M. Grossman, and J. C. Gee. Symmetric diffeomorphic image registration with cross-correlation: Evaluating automated labeling of elderly and neurodegenerative brain. *Med Image Anal*, 12(1):26–41, Feb 2008. [6](#)
- [12] Brian Avants, Jeffrey T Duda, Junghoon Kim, Hui Zhang, John Pluta, James C Gee, and John Whyte. Multivariate analysis of structural and diffusion imaging in traumatic brain injury. *Acad Radiol*, 15(11):1360–1375, Nov 2008. [4.2](#)
- [13] Brian Avants, Jeffrey T Duda, Junghoon Kim, Hui Zhang, John Pluta, James C Gee, and John Whyte. Multivariate analysis of structural and diffusion imaging in traumatic brain injury. *Acad Radiol*, 15(11):1360–1375, Nov 2008. [6](#)
- [14] Brian Avants, Alea Khan, Leo McCluskey, Lauren Elman, and Murray Grossman. Longitudinal cortical atrophy in amyotrophic lateral sclerosis with frontotemporal dementia. *Arch Neurol*, 66(1):138–139, Jan 2009. [4.2](#), [6](#)
- [15] Brian B Avants, Hallam Hurt, Joan M Giannetta, Charles L Epstein, David M Shera, Hengyi Rao, Jiongjong Wang, and James C Gee. Effects of heavy in utero cocaine exposure on adolescent caudate morphology. *Pediatr Neurol*, 37(4):275–279, Oct 2007. [4.2](#)
- [16] Brian B Avants, P. Thomas Schoenemann, and James C Gee. Lagrangian frame diffeomorphic image registration: Morphometric comparison of human and chimpanzee cortex. *Med Image Anal*, 10(3):397–412, Jun 2006. [2.7](#)
- [17] Sandhitsu R Das, Brian B Avants, Murray Grossman, and James C Gee. Registration based cortical thickness measurement. *Neuroimage*, 45(3):867–879, Apr 2009. [14](#), [3.3](#), [4.2](#), [6](#)
- [18] Jeffrey T Duda, Brian B Avants, Jane C Asmuth, Hui Zhang, Murray Grossman, and James C Gee. A fiber tractography based examination of neurodegeneration on language-network neuroanatomy. In *Workshop on Computational Diffusion MRI, Medical Image Computing and Computer-Assisted Intervention*, pages 191–198, Sep 2008. [6](#)
- [19] Jeffrey T Duda, Brian B Avants, Junghoon Kim, Hui Zhang, S Patel, John Whyte, and James C Gee. Multivariate analysis of thalamo-cortical connectivity loss in TBI. In *Computer Vision and Pattern Recognition*, pages 1–8, Los Alamitos, June 2008. IEEE Computer Society. [6](#)
- [20] M. Grossman, C. Anderson, A. Khan, B. Avants, L. Elman, and L. McCluskey. Impaired action knowledge in amyotrophic lateral sclerosis. *Neurology*, 71(18):1396–1401, 2008. [4.2](#)
- [21] M. Grossman, C. Anderson, A. Khan, B. Avants, L. Elman, and L. McCluskey. Impaired action knowledge in amyotrophic lateral sclerosis. *Neurology*, 71(18):1396–1401, Oct 2008. [6](#)
- [22] Junghoon Kim, Brian Avants, Sunil Patel, John Whyte, Branch H Coslett, John Pluta, John A Detre, and James C Gee. Structural consequences of diffuse traumatic brain injury: A large deformation tensor-based morphometry study. *Neuroimage*, 39(3):1014–1026, Feb 2008. [2.2](#), [2.7](#), [4.2](#)
- [23] Arno Klein, Jesper Andersson, Babak A Ardekani, John Ashburner, Brian Avants, Ming-Chang Chiang, Gary E Christensen, D. Louis Collins, James Gee, Pierre Hellier, Joo Hyun Song, Mark Jenkinson,

- Claude Lepage, Daniel Rueckert, Paul Thompson, Tom Vercauteren, Roger P Woods, J. John Mann, and Ramin V Parsey. Evaluation of 14 nonlinear deformation algorithms applied to human brain mri registration. *Neuroimage*, 46(3):786–802, Jul 2009. [1](#), [2](#), [4.2](#), [6](#)
- [24] Lauren Massimo, Chivon Powers, Peachie Moore, Luisa Vesely, Brian Avants, James Gee, David J Libon, and Murray Grossman. Neuroanatomy of apathy and disinhibition in frontotemporal lobar degeneration. *Dement Geriatr Cogn Disord*, 27(1):96–104, 2009. [4.2](#), [6](#)
- [25] J. Pluta, B. Avants, P. Yushkevich, S. Glynn, S. Awate, J. Detre, D. Mechanic, and J. C. Gee. Expectation matching for incomplete label driven semi-automated hippocampus segmentation in epilepsy. *Hippocampus*, page accepted, 2009. [1](#), [3.4](#), [4.2](#)
- [26] John Pluta, Brian B Avants, Simon Glynn, Suyash Awate, James C Gee, and John A Detre. Appearance and incomplete label matching for diffeomorphic template based hippocampus segmentation. *Hippocampus*, 19(6):565–571, Jun 2009. [6](#)
- [27] Tony J Simon, Zhongle Wu, Brian Avants, Hui Zhang, James C Gee, and Glenn T Stebbins. Atypical cortical connectivity and visuospatial cognitive impairments are related in children with chromosome 22q11.2 deletion syndrome. *Behav Brain Funct*, 4:25, 2008. [4.2](#)
- [28] Hui Sun, Brian B Avants, Alejandro F Frangi, Federico Sukno, James C Geel, and Paul A Yushkevich. Cardiac medial modeling and time-course heart wall thickness analysis. *Med Image Comput Comput Assist Interv Int Conf Med Image Comput Comput Assist Interv*, 11(Pt 2):766–773, 2008. [4.2](#)
- [29] Nicholas J Tustison, Brian B Avants, and James C Gee. Directly manipulated free-form deformation image registration. *IEEE Trans Image Process*, 18(3):624–635, Mar 2009. [4.2](#)
- [30] Paul A Yushkevich, Brian B Avants, John Pluta, Sandhitsu Das, David Minkoff, Dawn Mechanic-Hamilton, Simon Glynn, Stephen Pickup, Weixia Liu, James C Gee, Murray Grossman, and John A Detre. A high-resolution computational atlas of the human hippocampus from postmortem magnetic resonance imaging at 9.4 t. *Neuroimage*, 44(2):385–398, Jan 2009. [2.7](#), [4.2](#), [6](#)
- [31] Paul A Yushkevich, Brian B Avants, John Pluta, David Minkoff, John A Detre, Murray Grossman, and James C Gee. Shape-based alignment of hippocampal subfields: evaluation in postmortem mri. *Med Image Comput Comput Assist Interv Int Conf Med Image Comput Comput Assist Interv*, 11(Pt 1):510–517, 2008. [4.2](#)